User Guide

Run your Flex VRA monitors from your computer's keyboard



Get set up

1 Plug in the transmitter

Plug the USB transmitter into a USB-A port on the computer that runs your audiometry software. Flex PC finds it by itself.

2 Install Flex PC

Plug in the **silver USB stick** and run **FlexPC-Setup.exe** from it (no administrator password needed). It can add a desktop shortcut and start Flex PC with Windows.


3 Set your hotkeys

Click the gear to open **Settings**. Under **Hotkeys**, click a side's box and press a key (Ctrl, Shift or Alt allowed). Click **Save**.

No pairing needed. Every Flex PC USB comes pre-paired with your monitors. Just plug it in.

The Flex PC bar

Sits along the bottom of your screen. Drag it by the logo to move it.

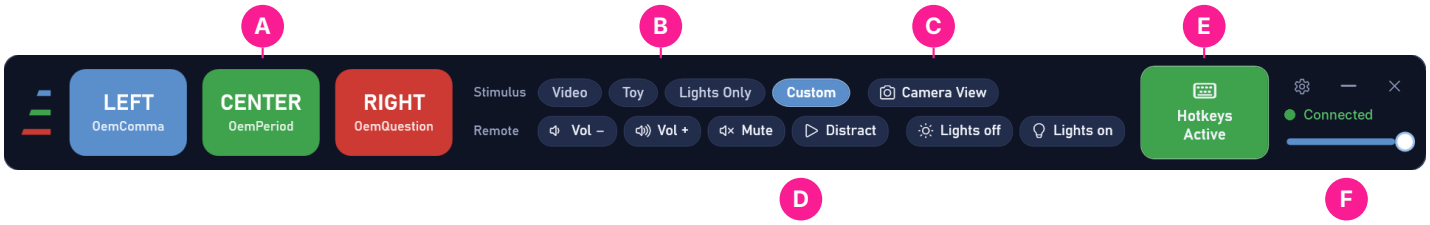

A **Left, Center, Right.** Hold the button or its hotkey to reinforce on that side.

B **Stimulus.** Video, Toy, Lights Only or Custom, as on the remote.

C **Camera View.** Shows the optional Flex Cam in the booth, and records the session.

D **Remote.** Volume, Mute, Distract, and the light bar on or off.

E **Hotkeys Active.** Click to pause the hotkeys while you type a report; click again to resume.

F **Settings and status.** Gear opens Settings. Green dot: transmitter connected. Slider: see-through.

> **Madsen Astera users:** Flex PC lets the buttons on your Astera audiometer's control panel run the VRA monitors (Tracking = Left, Reset = Center, Stim Lock = Right). **No pediatric module required.** Turn it on in **Settings › Astera control panel**.

During a test

- Keep your audiometry software in front as usual. The hotkeys work in any program, as long as the button shows **Hotkeys Active**.
- **Hold to activate** (the default): the side reinforces for as long as you hold the key, up to 30 seconds.
- **Timed**: one press reinforces for a set time, 1 to 30 seconds (4 unless changed). Choose this in **Settings › Activation**.

Settings

Activation	**Hold to activate** or **Timed** (seconds per press).
Hotkeys	A different key for each side. **Clear** removes one.
Astera control panel	Use the Astera panel's VRA buttons. Needs **Restart FlexPC as administrator**.
Transmitter	**Automatic**, or pick one. **Search again** looks for it.
Window	**Keep FlexPC on top of other programs** (on).

If something isn't working

The bar says Not connected	Unplug the transmitter and plug it back in, or open **Settings › Transmitter** and click **Search again**.
A hotkey does nothing	Check the button says **Hotkeys Active**, and that the key is set in **Settings › Hotkeys**.
The wrong monitor responds	Check the monitor's side in its own Setup (VRA page) matches the key you're pressing.

www.flexvra.com 512-777-0016 sales@flexvra.com Cedar Park, TX, USA